




Electric Mobile Stand

User Manual

<https://infraredredlight.com>

Adjustable Height

From 170cm (66.9") to 70cm (27.5")

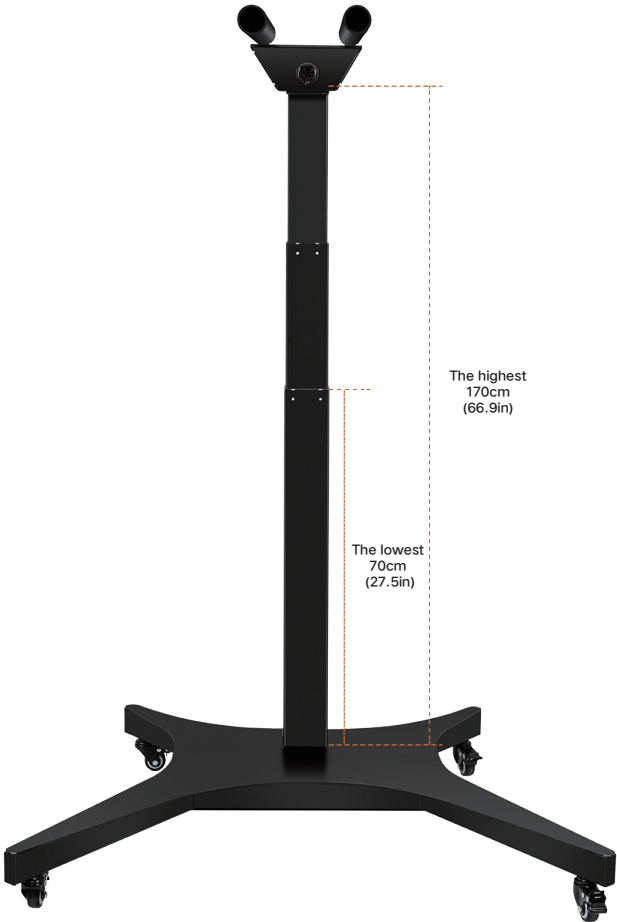

Component Parts




A) Square Post
(1)



B) Base
(1)



C) Universal Wheel
(4)

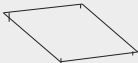

D) Base Plate
(4)



E) Remote Control
(1)



F) Power Cord
(1)



G) Adaptor
(1)



H) Control Box
(1)




I) Wire Controller
(1)



J) Fixing Rod For Lamp
(2)



K) Transverse Brace
(2)

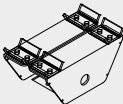

L) Lifting Rod Connection
(1)



M) Screws
M6*14(18)



N) Screws
M8*60(8)



O) 4mm Hexagon Wrench
(1)



P) 5mm Hexagon Wrench
(1)




Q) Cross Screws
M4*8(6)





R) Screws
M8*20(4)

1. Put the control system into the corresponding position according to the picture.

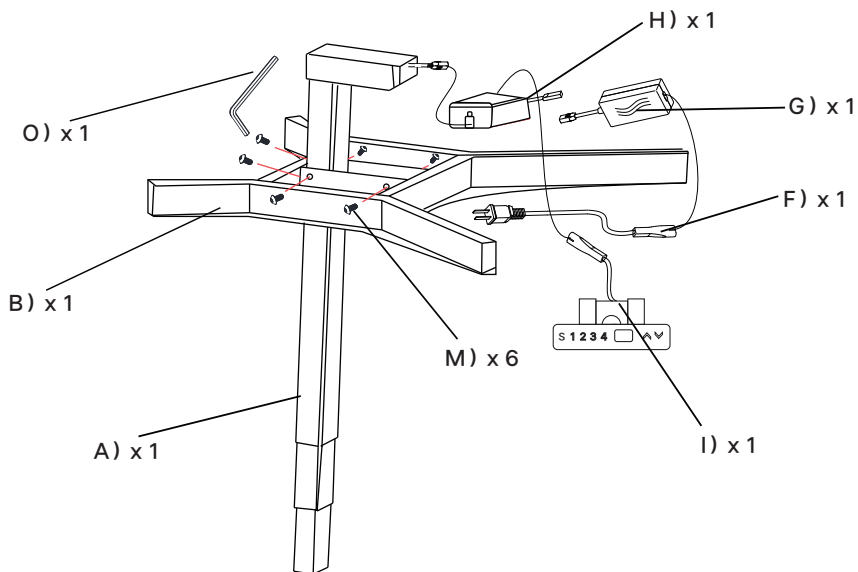


2. Fix cross screw M4*8 by cross screwdriver as picture.

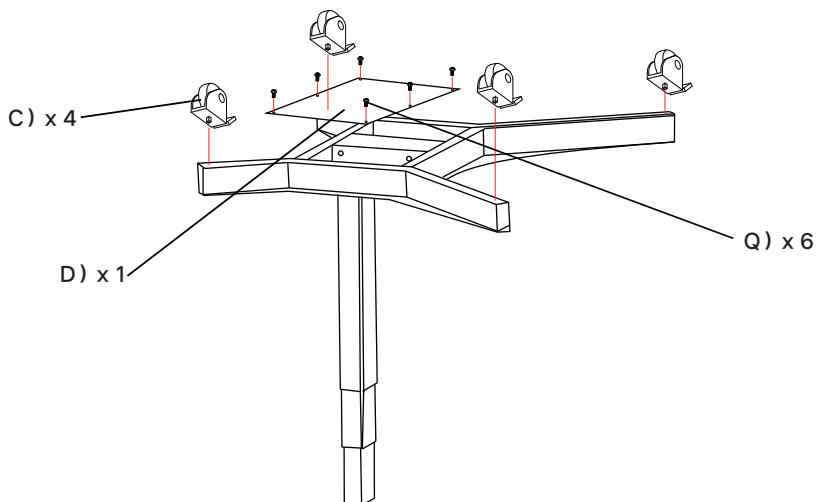

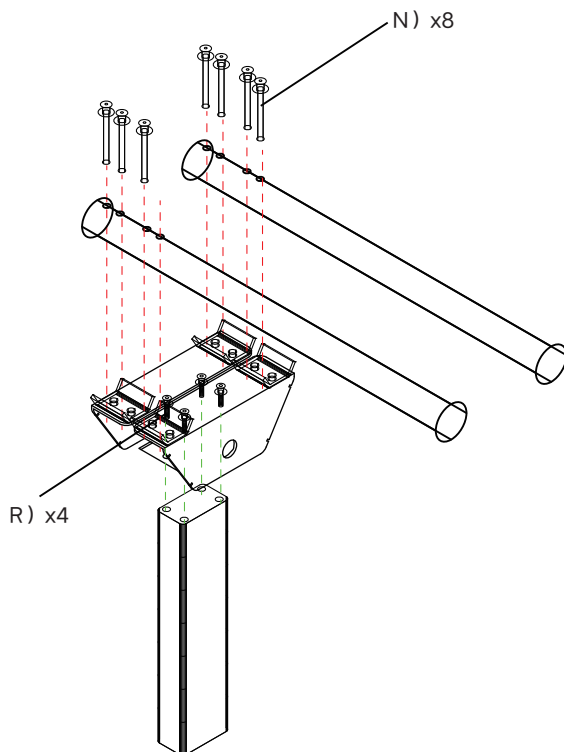


Controller Operation

Remote Control

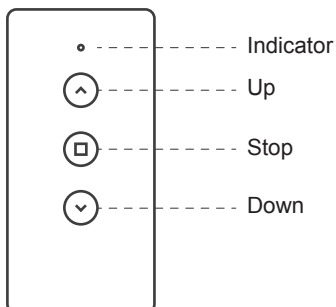

Wire Controller

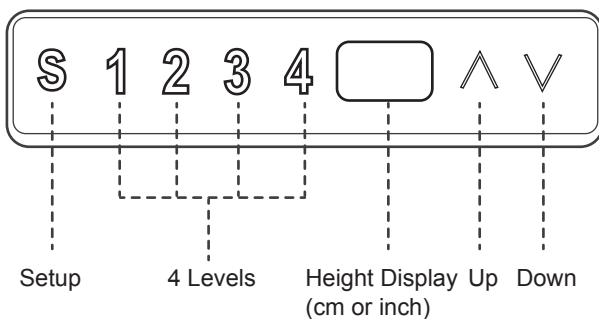

Switch the Height Unit

Press "S" button, if the digital shows "70-170", which means the height unit is centimeter (cm). You can switch the unit to an inch by lengthy pressing the "S" button until the digital shows "27.5-66.9".

Set the Height of Each Level

Adjust the height by pressing "^" and "V", then press the "S" button and one of "1-4" button to save the corresponding height of "1-4" gear. After setting, you can easily select four positions by pressing the "1-4" button as it applies memory function.

Adjust Height by Remote Controller

Press the button of "^" to rise, "V" to lower, "□" to stop rising and lowering.




Precautions

- 1) The product should be used properly.
- 2) Keep away from corrosive gas, liquid, and dusty environment.
- 3) Avoid continuously rising and lowering operation for over 1 minute, the display may show "HOT", it means the control system goes into a protection state, the remote and wire controller will be invalid correspondingly. Now, please turn off the power and wait about 20 minutes, then turn it back on and you can control it again.
- 4) After a prolonged use, it may happen that the height can not be risen, or the digital on the display shows "RST", then you need to reset it.



How to Reset?

Press the down button "V" to adjust the height to the lowest, then the height display of the wire controller shows a digital and the reset operation finishes. After that, you can use the controller to adjust the height normally.



Light Therapy Panels Mounting

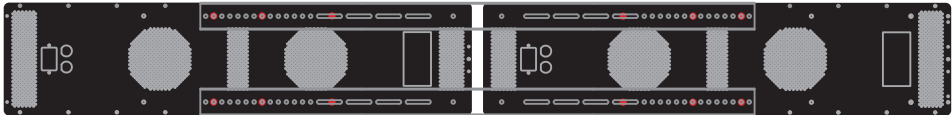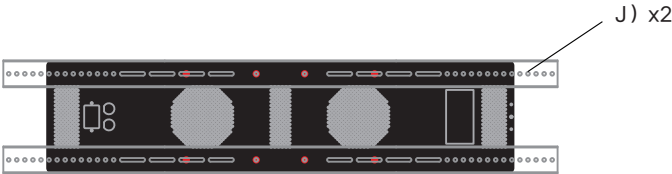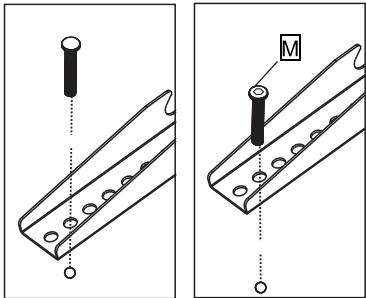


BP210

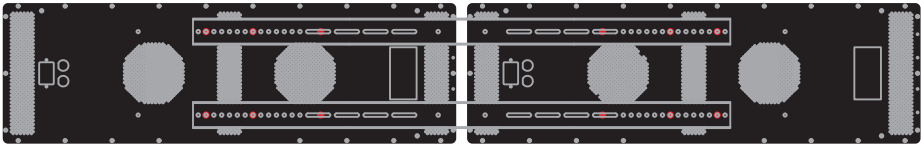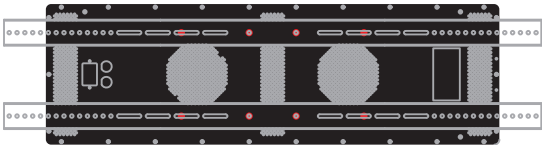



BP300

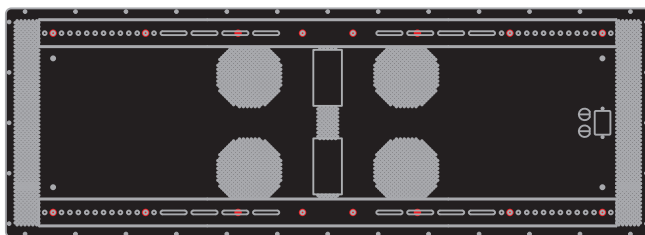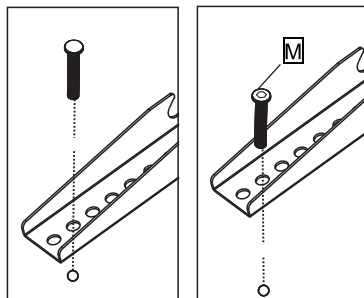

BP560

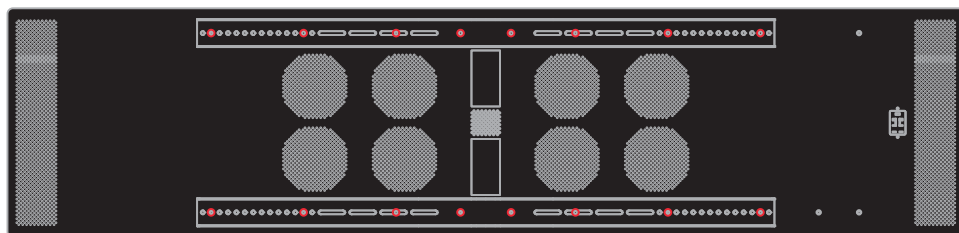

BP840

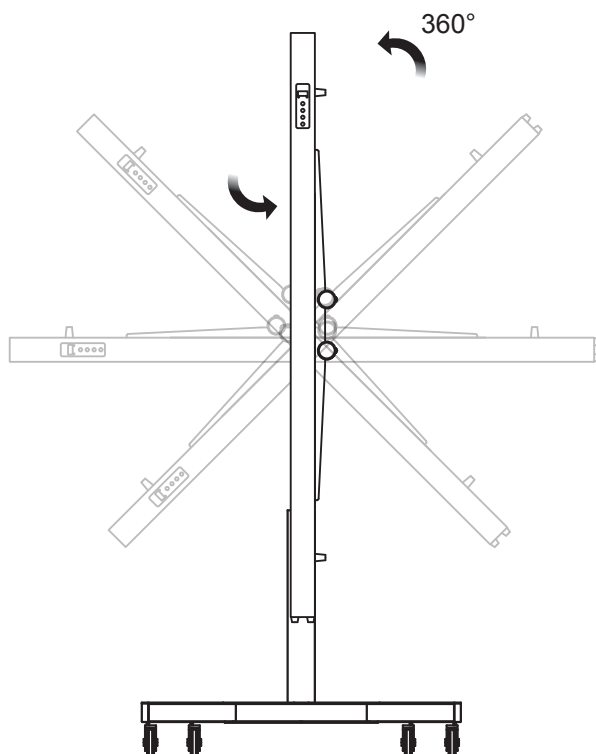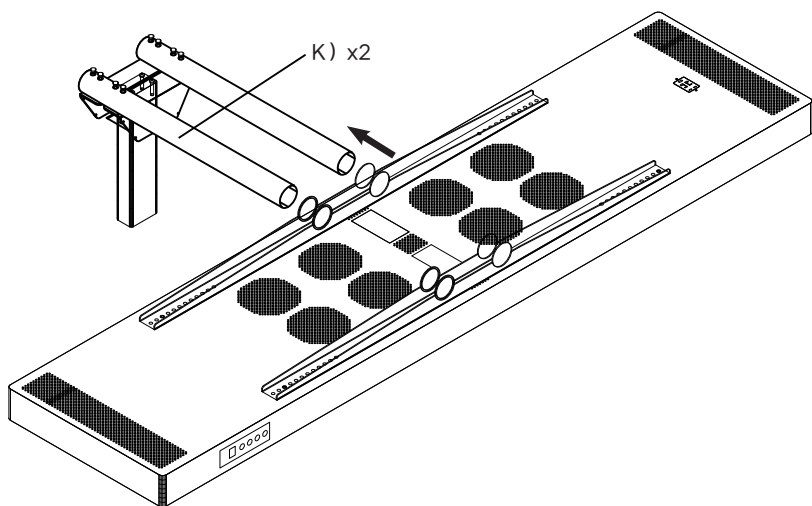

Thank You!